

# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

**RHW CLIN043**

Ref: T24/32297

<b>NAME OF DOCUMENT</b>	NAVA Nursing Management for Non-Invasive Ventilation
<b>TYPE OF DOCUMENT</b>	Clinical Business Rule
<b>DOCUMENT NUMBER</b>	RHW CLIN043
<b>DATE OF PUBLICATION</b>	20 May 2024
<b>RISK RATING</b>	Low
<b>REVIEW DATE</b>	May 2029
<b>FORMER REFERENCE(S)</b>	NAVA Nursing Management for Non-Invasive Mode
<b>EXECUTIVE SPONSOR</b>	S Bolisetty (Medical Co-Director Newborn Care Centre); S Wise (Nursing Co-Director Newborn Care Centre)
<b>AUTHOR</b>	E Deibe (ACNE), C Walter (CNE)
<b>SUMMARY</b>	To guide nurses on the care of neonates receiving non-invasive NAVA ventilation

# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

**RHW CLIN043**

This Clinical Business Rule is developed to guide safe clinical practice in Newborn Care Centre (NCC) at The Royal Hospital for Women. Individual patient circumstances may mean that practice diverges from this Clinical Business Rule. Using this document outside the Royal Hospital for Women or its reproduction in whole or part, is subject to acknowledgement that it is the property of NCC and is valid and applicable for use at the time of publication. NCC is not responsible for consequences that may develop from the use of this document outside NCC.

### 1. BACKGROUND

Neurally Adjusted Ventilatory Assist (NAVA) is a synchronised mode of ventilation intended for use in spontaneously breathing infants. NAVA can be used as both a Non-Invasive and Invasive mode via the Servo-n Getinge ventilator. NAVA modes require the placement of an Edi (electrical diaphragmatic index) catheter, which reads the diaphragmatic activity of the infant and adapts ventilation based on the measurement. This policy will guide nurses on the insertion of Edi catheter, interpretation of Edi monitoring and placement, understanding of NAVA levels and function, and the nursing care and management of infants on Non-Invasive NAVA ventilation.

### 2. RESPONSIBILITIES

Medical and Nursing Staff

### 3. PROCEDURE

#### 3.1 Equipment

- Getinge Servo-n Ventilator – set up and ready for use
- Fisher & Paykel snorkel
- Appropriate size hat and mask/prongs
- Edi catheter
- ENfit connector adaptor
- Leucoplast tape (thin)
- Barrier wipe
- Comfeel
- 1L bag of water for injection

#### 3.2 Clinical Practice

##### Commencing a patient on NAVA

1. Collect equipment.
2. Plug in Getinge ventilator to wall gas and power. Turn on (on/off switch is located at the back of the ventilator (Picture 1).



Picture 1


## NAVA Nursing Management for Non-Invasive Ventilation

RHW CLIN043

3. Commence pre-use check on ventilator by following instructions on screen.
4. Measure infants NEX score (Nare to tragus to xiphisternum) (Picture 2).

### NOTE:

- NEX measurement is NOT a deterrent of insertion depth for NG/OG tube.
- NEX is ONLY part of for the calculation of placement of Edi catheter.

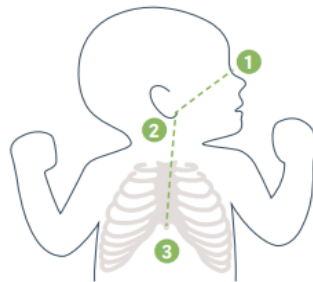

Picture 2

5. Click on Edi & NAVA on the main screen (Picture 3)
6. Select Edi catheter insertion calculation tool.
7. Calculate Edi catheter depth of insertion by entering all necessary information. (Picture 4)



Picture 3



Picture 4

8. Apply barrier wipe and Comfeel to cheek or chin to secure Edi catheter. Note that nasal insertion of Edi catheter is preferred.
9. Connect Edi catheter to Getinge Servo-n Edi module and cable.

# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

RHW CLIN043

10. Insert Edi catheter to pre-determined depth from calculation tool. Confirm placement on ventilator by using the positioning tool, assessing P wave formations (decreasing in size from top signal to bottom) (Picture 5).



Picture 5

11. Secure Edi catheter using brown tape. Note final depth of insertion.
12. Spike water bag with humidifier base line.
13. Turn on humidifier base once commencing patient use.
14. Ensure correct mode of ventilation has been chosen (NIV-NAVA/NIV-PC) and settings are correct as selected by the medical team.
15. Attach ENfit connector adaptor to allow oral syringe attachment to Edi catheter and confirm pH if feeding through.
16. Measure head circumference and choose appropriate size hat (Picture 6 & 7)
  - Hats should fit comfortably over the infant's head without leaving pressure marks or excessive gaping.
  - The hat should sit over the ears, round the back of the head and just above the eyebrows.



# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

RHW CLIN043

17. Select the correct size mask or prongs.

- The mask should provide an adequate seal over the nose/nostrils without causing pressure sores or excessive leak (Picture 8)
- The prongs should fill the nares without stretching the skin, but sit external to the nostril. The nose should not blanch when prongs are in position. (Picture 9)



Picture 8



Picture 9

18. Select appropriate size snorkel sizes.

19. Commence ventilation.

20. Place prongs/mask onto snorkel and secure to infant's face (Picture 10 & 11).

- Ensure straps are away from eyes and firm but not tight.
- You may need to alter the amount of foam under the snorkel to achieve an adequate seal and alleviate pressure on bridge of nose.



# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

**RHW CLIN043**

Picture 10

Picture 11

21. Document commencement of NIV NAVA ventilation on eRIC.

**NOTE:**

Ensure the Edi catheter is documented in eRIC.

- Click on *NICU Nursing* icon → *Procedures and Tubes* tab →select *NAVA* →click and follow prompts.
- Document the type of catheter (Gauge and length) as well as depth of insertion.

**Observations documented on eRIC**

- Continuous cardiorespiratory monitoring.
- Standard CPAP pressure cares and changes to be maintained 4<sup>th</sup> hourly.
- Blood gases PRN.

**Troubleshooting**

Problem/Alarm	Response
Edi catheter blockage	<p>Due to the addition of an electrical wire/probe, the Edi catheter has a smaller internal diameter lumen than the usual NG tube. This means that it can become blocked more easily. It is routine practice to have a second NG/OG tube inserted for the administration of feeds and medications.</p> <p>If the Edi catheter becomes blocked:</p> <ul style="list-style-type: none"> <li>• Remove the catheter, flush with 3-5ml of sterile water for injection and gently clean the tip of catheter with sterile water and gauze. Reinsert catheter to same depth.</li> <li>• If catheter remains blocked, remove, and replace with new Edi catheter.</li> </ul>
Excessive leak	<p>In non-invasive ventilation modes, a leak is to be expected. The ventilator will compensate for this and therefore leaks up to 90% may still guarantee effective ventilation.</p> <p>If “leakage too high” alarm:</p> <ul style="list-style-type: none"> <li>• Recheck fit of mask/prongs and snorkel to ensure snug (but not tight) fit. Resize or re-fit as needed.</li> <li>• Consider the use of chin strap to keep infants mouth closed- or use supported infant position to aim to keep mouth closed.</li> </ul>
No patient effort	<p>If machine alarming for No Patient Effort:</p> <ul style="list-style-type: none"> <li>• Ensure Edi is correctly positioned by checking depth of insertion at mouth/nare and Edi signal.</li> <li>• Check patient for respiratory effort to see if infant is breathing. If infant is breathing but Edi signal is</li> </ul>

# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

**RHW CLIN043**

	<p>saying infant is not breathing, there may be an issue with the electrodes.</p> <ul style="list-style-type: none"> <li>• Consider removing catheter, flushing with water and cleaning electrodes with gauze and sterile water before reinserting.</li> <li>• Back up rate should be initiated by machine if patient apnoeic.</li> <li>• Consider reducing the apnoea time if patient is not adequately supported in back up mode (discuss with medical team).</li> </ul>
<p>Gaseous distension</p>	<p>Because of the smaller internal diameter of catheter, difficulty in aspiration, as well as high flow rates with NIV NAVA, there may be air trapping in the gut leading to abdominal distention.</p> <ul style="list-style-type: none"> <li>• Place second NG/OG tube (if not already down) for aspiration of gas.</li> <li>• Monitor abdomen and immediately escalate any concerns including absence of bowel sounds, changes in abdominal colour, abdomen appearing ropey/veiny.</li> </ul>
<p>Volume Delivery Restricted Alarm</p>	<p>If you receive this alarm, it is usually because the infant is requiring a PIP (based on NAVA calculation) higher than what the ventilator will allow based on alarm settings. Medical staff should:</p> <ul style="list-style-type: none"> <li>• Enter alarm limits and reset PIP maximum alarm limit (as these functions as a pressure limiting alarm like PMAX). Note that the ventilator will cap the baby 5 below the set PIP max limit (i.e. if PIP alarm set at 25, maximum PIP that can be delivered is 20).</li> </ul>

### Pain and Comfort

- Ensure non-pharmacological interventions such as non-nutritive sucking, ISOC and dummy, comfort holding/kangaroo care and nesting/infant positioning are utilised for infant requiring NIV-NAVA.
- Prevent tug or drag from weight of tubing through positioning of infant, and supporting inspiratory and expiratory tubing.
- Attend 6 hourly CPAP cares to prevent pressure injuries and sores from developing.

### 3.1 Educational Notes

- Edi stands for electrical diaphragmatic index. This is measured in microvolts ( $\mu\text{V}$ ). When the respiratory centre in the brain is stimulated, it sends a signal via the phrenic nerve to the diaphragm muscle to stimulate contraction and subsequent respiration. When the muscle is stimulated it produces an electrical signal. The Edi catheter reads this electrical signal, displayed as Edi peak and minimum. Edi peak represents neural respiratory effort and is responsible for the size and duration of the breath. Edi min represents the spontaneous tonic (or resting) activity of

# Royal Hospital for Women (RHW)

## NEONATAL BUSINESS RULE

### NAVA Nursing Management for Non-Invasive Ventilation

RHW CLIN043

the diaphragm in between inspirations, which prevents de-recruitment of the alveoli during expiration.

- Edi Peak is 5-15.
  - If Edi peak is >15, this indicates the diaphragm muscle is working hard and that the babies breathing is being under-supported. This can be seen clinically through increased work of breathing, increased oxygen requirements, more frequent desaturations or bradycardias or increasing respiratory acidosis on a gas.
  - If Edi peak is <5 there is a weak signal from brain to diaphragm and is a sign that the baby is being over-supported. This can also present as an increase in frequency of apnoeic episodes due to hypocapnia.
  - In an infant requiring higher percentage of back-up and having increasing periods of 'no patient effort' it is important to assess what the Edi peak trends have been as this could be a sign of the baby being under or over supported .
- Edi min
  - Edi min should be targeted <3.
  - If Edi min is consistently greater or equal to 3, consider increasing PEEP to reduce the tonic activity of the diaphragm and to maintain Functional Residual Capacity.
- Edi Trigger
  - Edi trigger is the minimum increase in electrical activity that triggers the ventilator.
  - Usually set at 0.5 microvolts and triggers NAVA to assist with the breath.
  - If the Edi trigger is set too low, the ventilator responds to small Edi signals and converts them into small breaths. That prevents neonates from going into backup ventilation and may result in under-ventilation.
- NAVA Level
  - Determines the positive inspiratory support the baby receives based on the Edi peak and min values. PIP level is determined based on the following calculation:

$$(Edi\ Peak - Edi\ Min) \times NAVA + PEEP$$

  - For example, if NAVA level is set at 1.0, Edi Peak is recorded as 15, Edi min was 2 and PEEP is set at 6 cm then the ventilator calculates  $(15-2) \times 1 + 6 = 19$ . This means on the next breath, the ventilator will support infant to a PIP of 19. The neonate is also responsible for determining the inspiratory and expiratory times and respiratory rate whilst actively breathing.
  - If the infant does not make respiratory effort, however, the calculated PIP based on the NAVA level is disregarded and the infant will automatically receive the pre-set back up settings.
- NAVA Level 0
  - In some patients, the medical team will choose to place baby on NIV NAVA with a NAVA level of 0. This means that no additional PIP support (above PEEP) is given when the infant is actively breathing – providing CPAP pressures. Using the above example  $(15-2) \times 0 + 6 = 6$ .
  - The benefit of nursing an infant on NAVA level of 0 is that if patient becomes apnoeic and has no respiratory effort, back up mode will still be triggered. Therefore, NIV NAVA with a NAVA level of 0 is the same as giving CPAP with back up.
- Apnoea Time
  - Determines the amount of time the neonate can be apnoeic before the ventilator switches into backup mode. Apnoea is typically defined as no respiratory effort for 20 seconds, however this typically results in clinical symptoms including bradycardia and desaturation.
  - Setting a shorter apnoea time allows the machine to revert into back up, preventing clinical signs of apnoea.



# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

RHW CLIN043

- Apnoea time is generally set at 2 seconds i.e. after 2 seconds of no respiratory effort, the ventilator will switch into back-up mode.
- A neonate who remains apnoeic will continue to be ventilated in backup mode until respiratory effort returns.
- NIV PC
  - NIV Pressure control. Pre-set inspiratory pressure is delivered at a pre-set respiratory rate.
- Other considerations
  - Flexitrunk nasal tubing (clear, rigid part of the snorkel) should be the shortest length possible; not extending over the infant's head and parallel to face, using foam to adjust if needed.
  - Internal expiratory block of ventilator only needs to be changed at the end of treatment when it has been used on an infected patient (e.g. MRSA, serratia, RSV). At all other times it can be reused if the additional disposable filter has been in situ for the duration of the treatment.
  - The disposable filter attached to white expiratory line should be changed 48 hourly during use
  - Ventilator circuit should be changed every 14 days.
- Edi Tips
  - Ensuring the Edi catheter is in the correct place is essential to ensuring that the patient is receiving effective ventilation, given that PIP is determined by catheter readings.
  - Edi signal should be checked every hour at a minimum, and any changes to signal should be rectified promptly.
  - Looking at the Edi placement screen, the waves on the top should have pronounced P & QRS complex. As you move down the screen dampened or absent P waves and dampened QRS complex appears indicating ideal placement. On smaller babies, it is not uncommon that P waves will still be pronounced in the bottom waveform due to span of electrodes and cardiac interference. Therefore the aim should be to see some decreasing size of P waves from top to bottom, even if absence of P wave in bottom waveform can't be achieved.
  - The catheter insertion positioning tool will indicate if movement of the Edi catheter needs to occur (Picture 12 and 13).



Picture 12



Picture 13

# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

**RHW CLIN043**

- Edi catheters can be used for 4 weeks. The site should be rotated weekly, and the catheter cleaned with sterile water and dried with gauze prior to reinsertion. More frequent changes may be required if feeding through catheter.
- Edi catheters can be placed orally or nasally. Be aware of the risk of pressure injuries that can occur, in particular, on very small infants. Despite preference for nasal insertion (due to less movement of the catheter and therefore more reliable readings), an orally placed Edi may be necessary for small infants.

### 3.4 Abbreviations

NAVA	Neurally Adjusted Ventilatory Assist	PRN	Pro re neta (as required)
Edi	Electrical Diaphragmatic Index	eRIC	Electronic Record for Intensive Care
NEX	Nare, Ear, Xiphisternum	PIP	Positive Inspiratory Pressure
NG	Nasogastric	ISOC	Immuno-supportive Oral Care
OG	Orogastric	CPAP	Continuous Positive Airway Pressure
PMAX	Pressure Maximum	µV	Micro Volt
NIV-NAVA	Non- Invasive Neurally Adjusted Ventilatory Assist	PEEP	Positive End Expiratory Pressure
NIV-PC	Non- Invasive Pressure Control	MRSA	Methicillin-resistant Staphylococcus Aureus
pH	Power of Hydrogen	RSV	Respiratory Syncytial Virus
NICU	Neonatal Intensive Care Unit		

### 3.5 References

1. Maquet Geringe Group (2015) Servo-n Self-Guided Education Presentations. Maquet. Rastatt (Germany)
2. Maquet Geringe Group (2015) Neurally Adjusted Ventilatory Assist (NAVA) - Synchrony redefined. [online] available from [www.maquet.com](http://www.maquet.com) (Accessed on 3/12/23)
3. Maquet Geringe Group (2013) Ventilation Servo-I for Neonates. Synchrony for those who need it most. Maquet. Solna (Sweden)
4. Servo-n Hands on guide. 2023. <https://getinge.training/wp-content/uploads/2023/02/MCV00096728-REVA-Servo-n-2.1-Training-Hands-On-Guide-for-US.pdf> Accessed 24/03/2024.

## 4. RELATED BUSINESS RULES AND POLICY DOCUMENTS

- NAVA Nursing Care for Invasive Ventilation
- RHW NCC Nursing- Geringe Servo-n set up
- RHW NCC Medical- NAVA (Neurally Adjusted Ventilatory Assist) Clinical Guidelines

## 5. CULTURAL SUPPORT

- When clinical risks are identified for an Aboriginal family, they may require additional supports. This may include Aboriginal health professionals such as Aboriginal liaison officers, health workers or other culturally specific services.

# Royal Hospital for Women (RHW) NEONATAL BUSINESS RULE

## NAVA Nursing Management for Non-Invasive Ventilation

**RHW CLIN043**

- For a Culturally and Linguistically Diverse CALD family, notify the nominated cross-cultural health worker during Monday to Friday business hours.
- If the family is from a non-English speaking background, call the interpreter service: NSW Ministry of Health Policy Directive PD2017\_044-Interpreters Standard Procedures for Working with Health Care Interpreters.

### 6. IMPLEMENTATION PLAN

This revised CBR will be distributed to all medical, nursing and midwifery staff via @health email. The CBR will be discussed at ward meetings, education and patient quality and safety meetings. Education will occur through in-services, open forum and local ward implementation strategies to address changes to practice. The staff are asked to respond to an email or sign an audit sheet in their clinical area to acknowledge they have read and understood the revised CBR. The CBR will be uploaded to the CBR tab on the intranet and staff are informed how to access.

### 7. RISK RATING

- Low

### 8. NATIONAL STANDARDS

- Standard 3 Preventing and Controlling Infections
- Standard 5 Comprehensive Care
- Standard 6 Communicating for Safety
- Standard 8 Recognising and Responding to Acute Deterioration

### 9. REVISION AND APPROVAL HISTORY

Date	Revision No.	Author and Approval
04/10/2016	1	A Ottoway (ACNE/CNS)
20/03/2024	2	E Deibe (ACNE/CNS), C Walter (CNE) Approval NCC CBR Committee
23.4.24		Endorsed at BRGC